

Zenex Multi & Plus Kit


Zenex Prosthetic Kit



IZ-CAT-05 Rev.00 (06/20)

Z E N E X MULTI&PLUS K I T

012 Zenex Prosthetic Kit



IZEN Implant, innovative technology in pursuit of **bright smile**


Established in 2016 by developers and researchers who have spent more than 20 years on research and development in dental implants, IZEN Implant have working on developing innovative products to provide the best result to patients. We will continue to grow into a company that meets customers' needs with safe and excellent products made of trust and continuous research and development.

IZEN Implant provides more convenient, stable, and optimized implant system by building product line-up with high-end technology of Zenex Multi and Zenex Plus.

Certification



History & Global Network

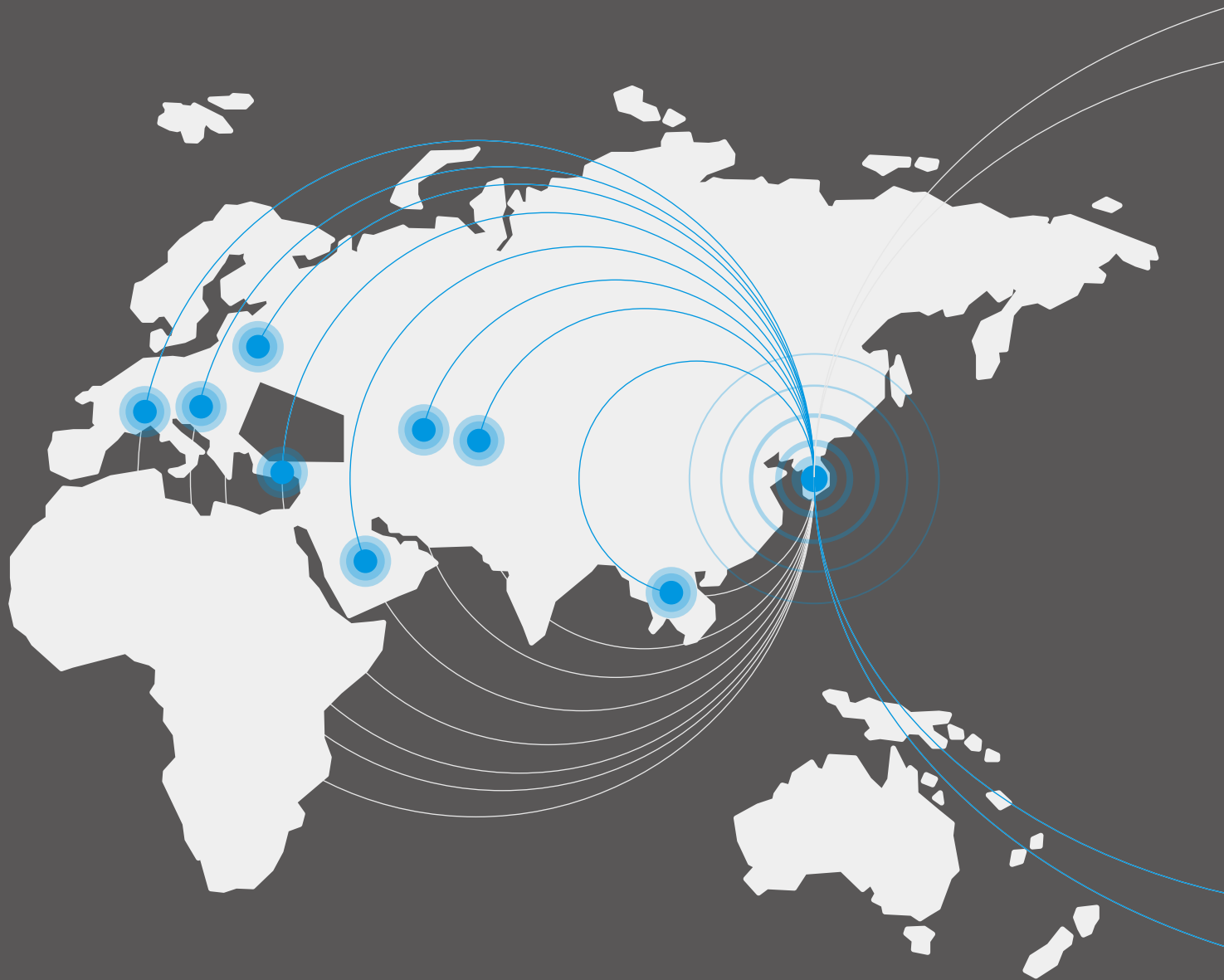

2016. 03

Foundation of IZEN Implant

2016. 06

Acquisition of manufacturing licenses for medical devices
Report medical device items
(Screw driver for dental implant surgery, No. 16-760)

2016. 12

Patent registration
(Driver holder and driver apparatus for dental engineering, No. 10-1684959)

2017. 03

Patent registration
(Driver holder for dental apparatus, No. 10-1714533, No. 10-1714541)

2018. 08

Venture Business Certificate
(No. 20180110278)

2018. 09

Foundation of Corporate
Affiliated Research Institute (No. 2018114470)

2019. 06

Acquisition of small businesses
technology development project
(Development of Dental implant abutments and scan body for

2019. 11

Acquisition of ISO13485:2016
(GIS-4058-MD)



2017. 09

Registration of medical devices
(Apparatus for dental implant
surgery, No. 17-1033)

2018. 03

Certification of medical
device items compliance with
manufacturing and quality
control standards (KGMP)

2018. 04

Medical Device Item Certification
(Dental implant abutments, No.
18-4268)

2018. 07

Acquisition of ISO 9001 : 2015
(K-GM-4699)

2019. 12

Patent registration
(Abutment and dental implant,
No. 10-2062575)

Kit Contents

Zenex
Prosthetic Kit
012



1.2 Hex
Driver
014



1.2 Hex Machine
Screw Driver
014



Abutment
Remover
015



Finger Driver
for Driver Shaft
015



Ball Abutment
Driver
016



Multi Straight
Abutment Driver
016



Torque
Wrench
017



Kerator
Magic Tool
017



Kerator
Torque Tip
018



IZEN
I M P L A N T

009
ZENEX
KIT

IZEN


I M P L A N T

ZENEX PROSTHETIC KIT

012 ZENEX
KIT

Zenex Prosthetic Kit

Composition on lower kit tray

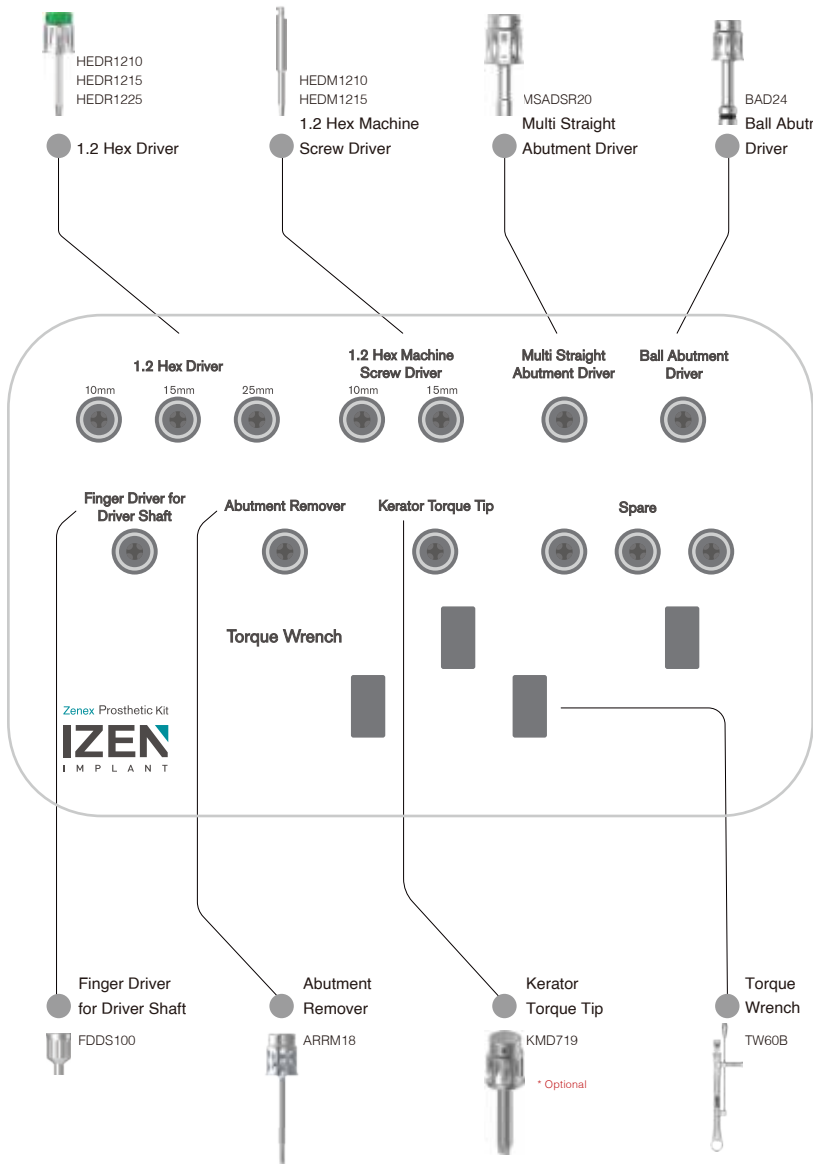


Zenex Prosthetic Kit Instruments

1.2 Hex Driver

- Driver for connecting with Torque Wrench
- Used for connecting to cover screw and surgical instrument with 1.2 Hex
- Follow recommended tightening torque calculation value
- Be sure to check transformation of Hex tip before use

L (mm)	10.0	15.0	25.0
	HEDR1210	HEDR1215	HEDR1225

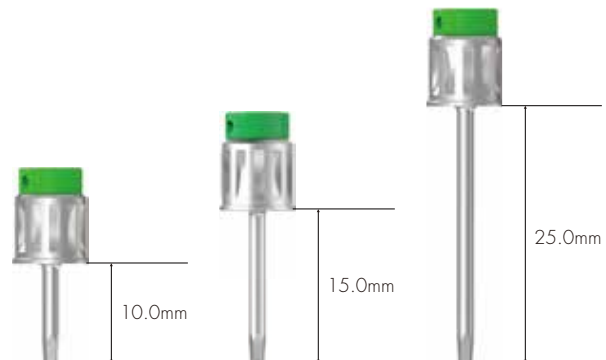


014
ZENEX
KIT

1.2 Hex Machine Screw Driver

- Driver for implant motor
- Implement holding instrument
- Available to use by hand when connected to Finger Driver for Driver Shaft

L (mm)	10.0	15.0
	HEDM1210	HEDM1215

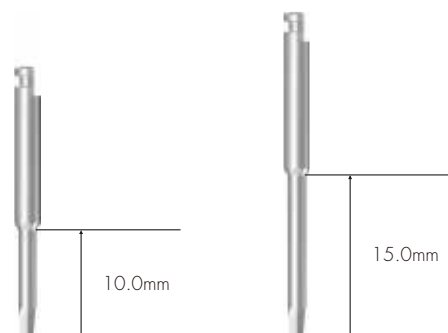

Abutment Remover

- Use for removing abutment when it is impossible to remove abutment by hands due to strong connection force.
 - Remove Abutment by keep tightening the remover clockwise after removing abutment screw
- ※ R-System use only



ARRM18

Finger Driver for Driver Shaft

- Available to rotate with hand by connecting to every surgery instrument for motor



FDDS100

Zenex Prosthetic Kit Instruments

Ball Abutment Driver

- Torque driver for Ball Abutment



BAD24

016 ZENEX
KIT

Multi Straight Abutment Driver

- Torque driver for Multi Straight Abutment



MSADSR20

Torque Wrench

- Used for placing fixture and tightening implant abutments and screws
- Confirm torque calculation value by pulling the bar to be aligned with the marked line
- Recommended to use by pressing head to prevent detachment of Torque Wrench and surgical instruments



TW60B

Kerator Magic Tool

- Use to connect and remove replacement male on denture cap
- As separated into two pieces, it is easy for maintenance (Optional)



KMH002

Kerator Torque Tip

- Torque driver for Kerator Abutment (Optional)

Round Type



KMD719

IZEN

I M P L A N T

019
ZENEX
KIT