



ZENEX Implant System

CONTENTS

01	Zenex Mini Submerged System	
	Zenex Mini Submerged Fixture	006
	Cover Screw	006
	Prosthetic Flow Diagrams	007
	Healing Abutment	008
	Temporary Abutment	008
	Simplex Abutment	009
	Duplex Abutment	010
	Angled Abutment	011
	Pick-up Impression Coping	012
	Transfer Impression Coping	012
	Fixture Lab Analog	012
02	Zenex Mini Post System	
	Zenex Mini Post Fixture	014
	Prosthetic Flow Diagrams	016
	Mini Impression Coping	017
	Temporary Cap	017
	Post Lab Analog	017
03	Zenex Mini Ball System	
	Zenex Mini Ball Fixture	020
	Prosthetic Flow Diagrams	022
	Ball Lab Analog	023
	Ball Cap O-ring Set	023

01


ZENEX Mini Submerged System

ZENEX Mini Submerged Fixture

Screw & Cement Retained Restoration

ZENEX Mini Submerged (Two Piece)

Diameter Ø3.2

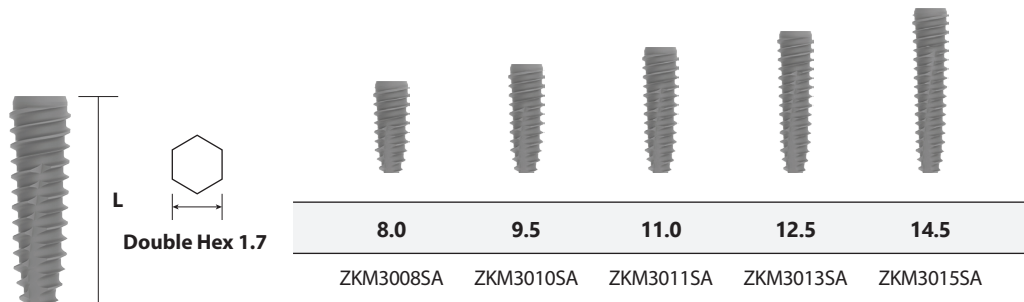

• Recommendable torque for Placement : Less than 40Ncm

+FDA approval is under processing

Cover Screw




ZKMCS14

- Use 1.25 Hex driver
- Torque to tighten : 8 ~ 10Ncm
- Screw - M1.4


Prosthetics Flow Diagrams for ZENEX Mini Submerged System

Screw & Cement Retained Restoration

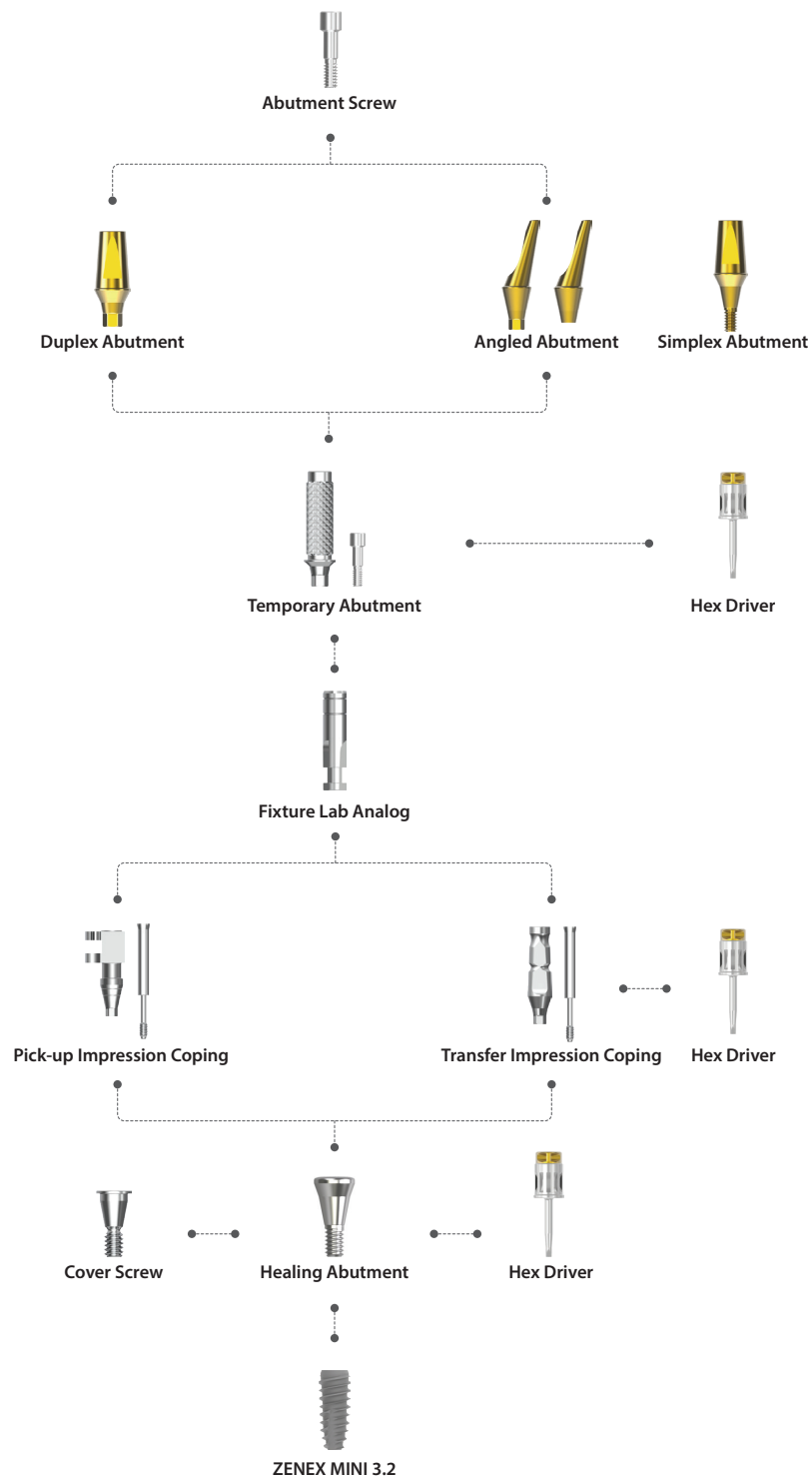


Healing Abutment

Ø3.0

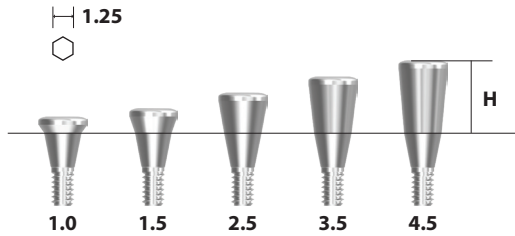


Ø3.5

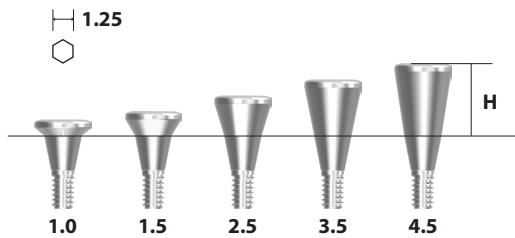

L \ D	Ø3.0	Ø3.5
1.0	ZKMHA3010	ZKMHA3510
1.5	ZKMHA3015	ZKMHA3515
2.5	ZKMHA3025	ZKMHA3525
3.5	ZKMHA3035	ZKMHA3535
4.5	ZKMHA3045	ZKMHA3545

- Use 1.25 Hex driver
- Torque to tighten : 5 ~ 8Ncm
- Screw - M1.4


Temporary Abutment

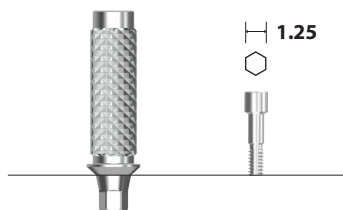

P/H	G/H \ D	Ø3.5
1.0	10.0	ZKMTA35H

- Use to make the temporary prosthesis
- The structure is easy to customize
- Use 1.25 Hex Driver
- Torque to tighten : 20Ncm
- Abutment Hex 1.7 / Screw - M1.4


Simplex Abutment

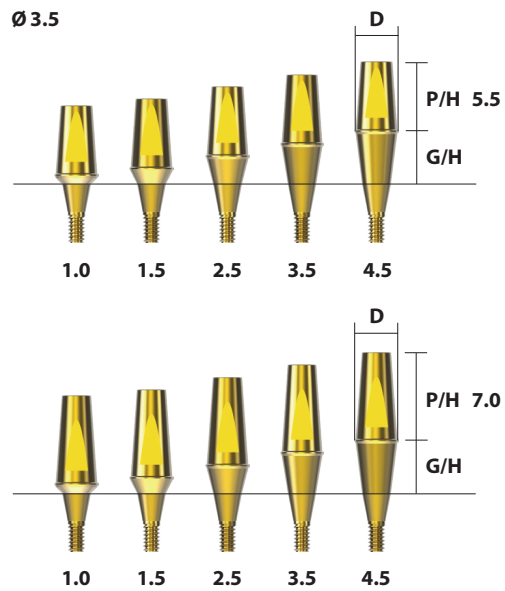

P/H	G/H \ D	Ø3.5
5.5	1.0	ZKMSA351055
	1.5	ZKMSA351555
	2.5	ZKMSA352555
	3.5	ZKMSA353555
	4.5	ZKMSA354555
7.0	1.0	ZKMSA351070
	1.5	ZKMSA351570
	2.5	ZKMSA352570
	3.5	ZKMSA353570
	4.5	ZKMSA354570

- Designed for general cement type prosthesis
- One body structure (abutment with screw)
- 11° taper structure for stable connection
- Colored in gold for esthetic satisfaction
- Profile for anti-rotation of prosthesis
- Use 1.25 Hex driver
- Torque to tighten : 30Ncm
- Screw -M1.4


Duplex Abutment

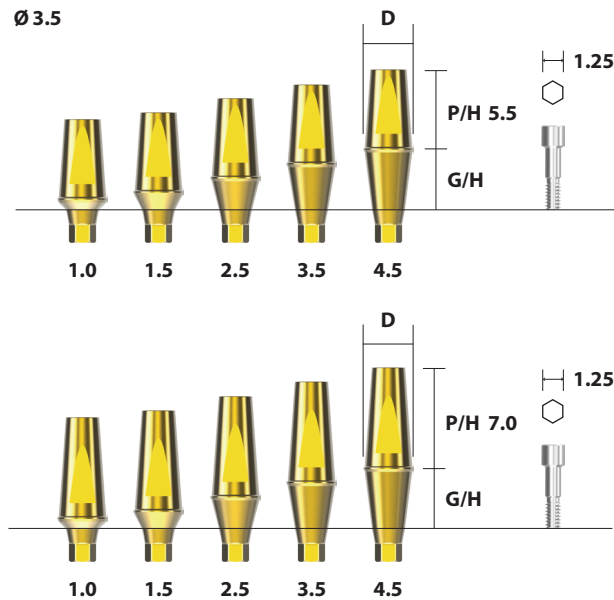

P/H	G/H \ D	Ø3.5
5.5	1.0	ZKMDA351055H
	1.5	ZKMDA351555H
	2.5	ZKMDA352555H
	3.5	ZKMDA353555H
	4.5	ZKMDA354555H
7.0	1.0	ZKMDA351070H
	1.5	ZKMDA351570H
	2.5	ZKMDA352570H
	3.5	ZKMDA353570H
	4.5	ZKMDA354570H
Screw		ZKMAS14

- Designed for general cement type prosthesis
- One body structure (abutment with screw)
- 11° taper structure for stable connection
- Colored in gold for esthetic satisfaction
- Profile for anti-rotation of prosthesis
- Use 1.25 Hex driver
- Torque to tighten : 20Ncm
- Abutment Hex 1.7 / Screw - M1.4

Angled Abutment

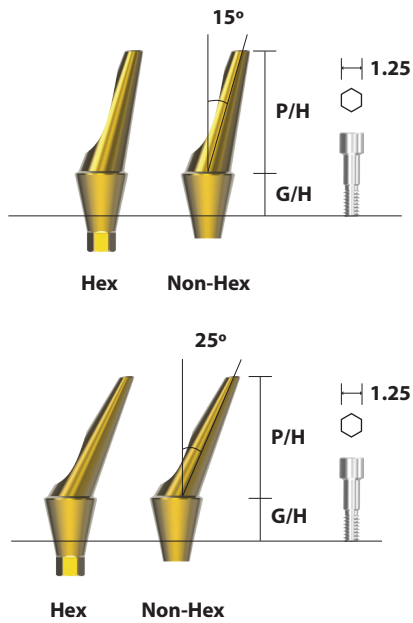

D			Ø3.5	
Angled	P/H	G/H \ Type	Hex	Non-Hex
15°	8.0	2.0	ZKMAA35215H	ZKMAA35215N
		3.0	ZKMAA35315H	ZKMAA35315N
		4.0	ZKMAA35415H	ZKMAA35415N
25°	8.0	2.0	ZKMAA35225H	ZKMAA35225N
		3.0	ZKMAA35325H	ZKMAA35325N
		4.0	ZKMAA35425H	ZKMAA35425N
Screw			ZKMAS14	

- It is applied to adjust a path for prosthesis with its 15° & 25° axial angle
- 11° tapered structure for stable connection
- Colored in gold for esthetic satisfaction
- Use 1.25 Hex driver
- Torque to tighten : 20Ncm
- Abutment Hex 1.7 / Screw - M1.4

Pick-up Impression Coping



L \ D	Ø3.5
13.3	ZKMIP35H
Guide Pin	ZKMIGS
	ZKMIGL

+ Option : Short Guide pin / ZKMIGS

- Pick-up type impression coping with the custom tray
- The item is designed to minimize ill-fitting application in the impression material
- Asymmetric design is applied to minimize adjacent interference
- Height variations (Long & Short) of Guide Pin for convenience
- Use 1.25 Hex Driver
- Abutment Hex 1.7 / Screw -M1.4

Transfer Impression Coping



L \ D	Ø3.5
14.3	ZKMIT35H
Guide Pin	ZKMIGS
	ZKMIGL

+ Option : Short Guide pin / ZKMIGS

- Transfer type impression coping with the custom tray
- Curved contour for increasing insertion ability
- Height variations (Long & Short) of Guide Pin for convenience
- Use 1.25 Hex Driver
- Abutment Hex 1.7 / Screw -M1.4

Fixture Lab Analog



ZKMFLA30

- The item which makes on the model of the actual fixture

02

ZENEX Mini
Post System


ZENEX Mini Post Fixture

Screw & Cement Retained Restoration


ZENEX Mini (Post)

Diameter Ø2.8

G/H = 2.5

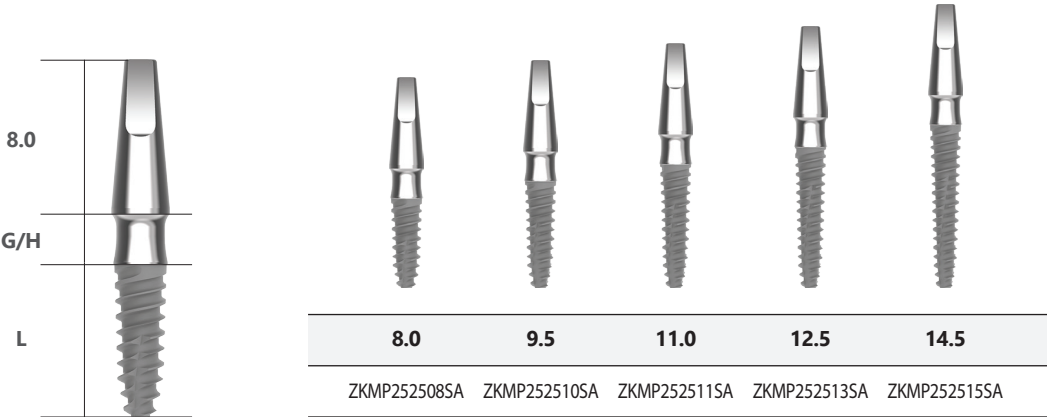

G/H = 4.0

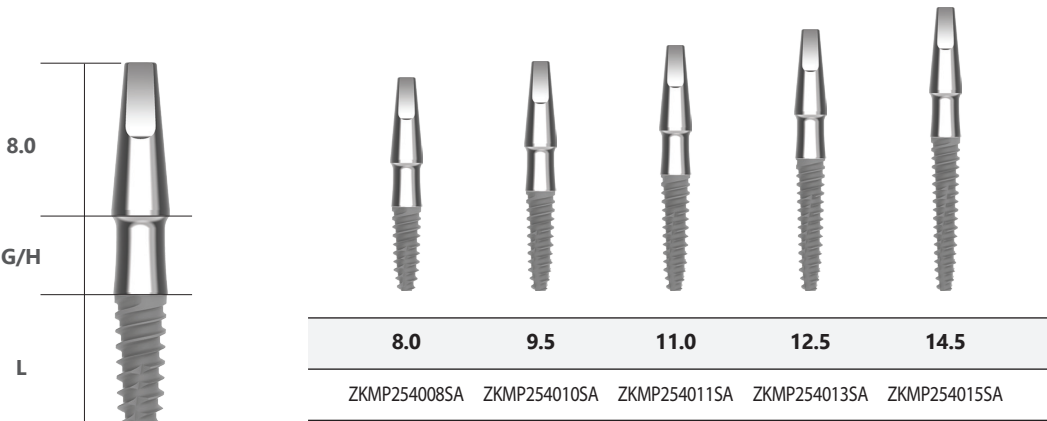


• Recommendable torque for Placement : Less than 40Ncm
 + FDA approval is under processing

ZENEX Mini (Post)


Diameter Ø3.0

G/H = 2.5

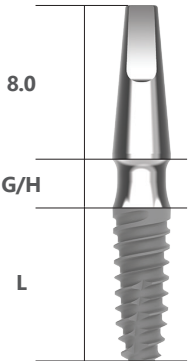

				
8.0	9.5	11.0	12.5	14.5
ZKMP302508SA	ZKMP302510SA	ZKMP302511SA	ZKMP302513SA	ZKMP302515SA

G/H = 4.0

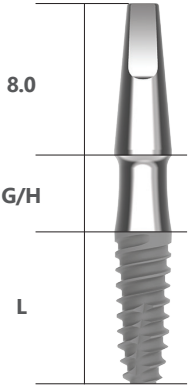


				
8.0	9.5	11.0	12.5	14.5
ZKMP304008SA	ZKMP304010SA	ZKMP304011SA	ZKMP304013SA	ZKMP304015SA

• Recommendable torque for Placement : Less than 40Ncm
+ FDA approval is under processing

Prosthetics Flow Diagrams for ZENEX Mini Post System

Cement Retained Restoration / Overdenture Retained Restoration

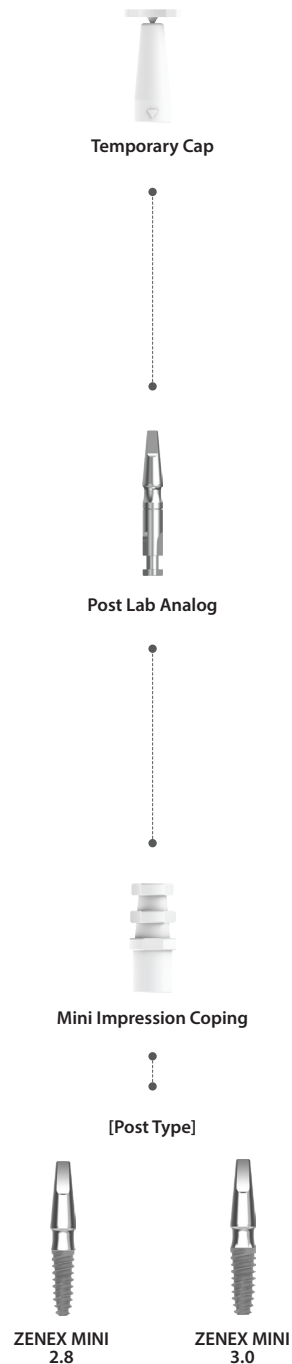

Mini Impression Coping



ZKMIOP

- Used in taking impression from the Zenex Mini fixture (Post Type)
- Screw -M1.4

Temporary Cap



ZKMTC

- Use it to make the temporary prosthesis

Post Lab Analog



ZKMPLA

- The item which makes on the model of the actual Zenex Mini fixture (Post Type).

03


**ZENEX Mini
Ball System**

ZENEX Mini Ball Fixture


ZENEX Mini (Ball)

Diameter Ø2.8

G/H = 2.0

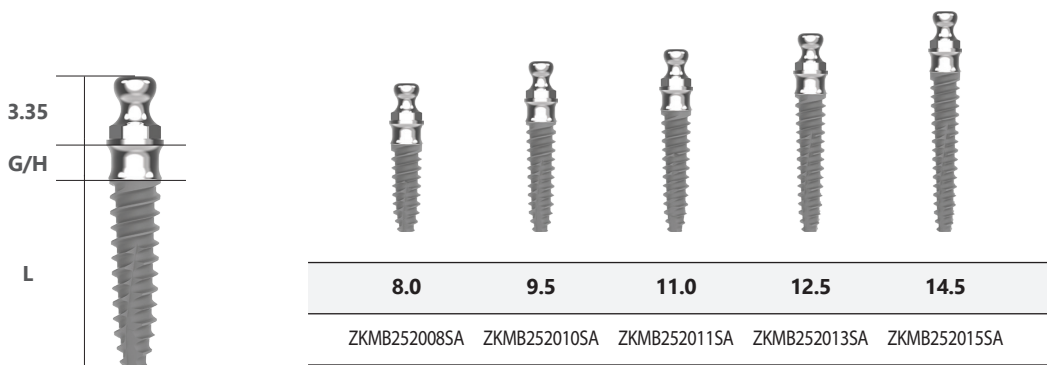

G/H = 4.0

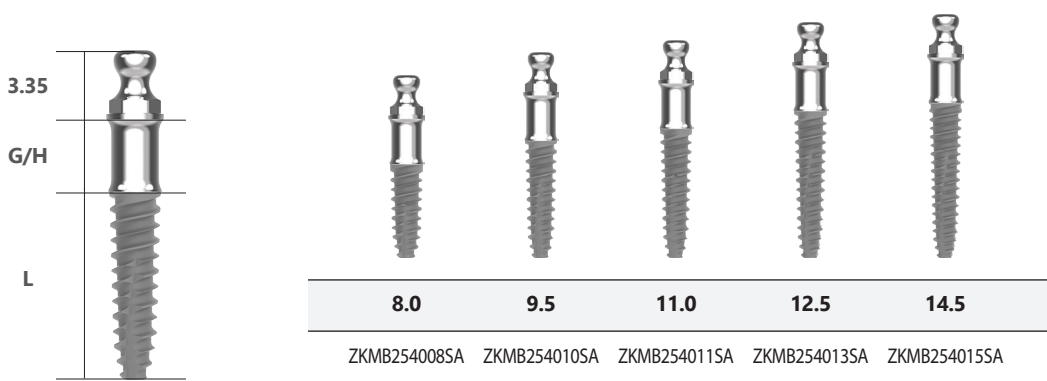


•Recommendable torque for Placement : Less than 40Ncm
+ FDA approval is under processing


ZENEX Mini (Ball)

Diameter Ø3.0

G/H = 2.0

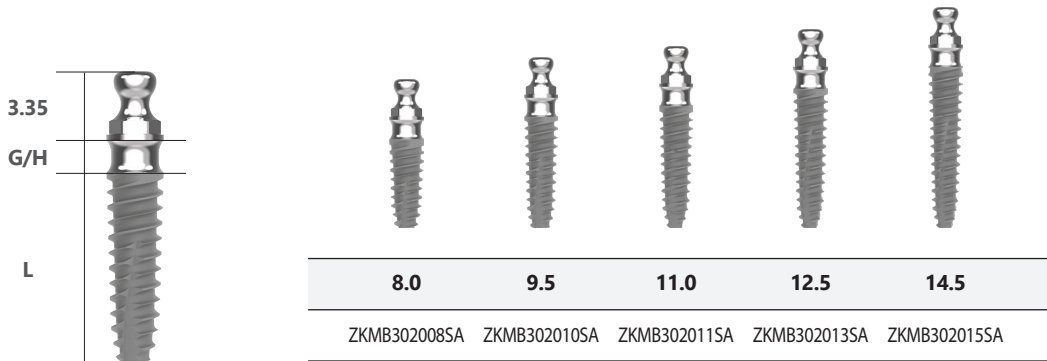


G/H = 4.0

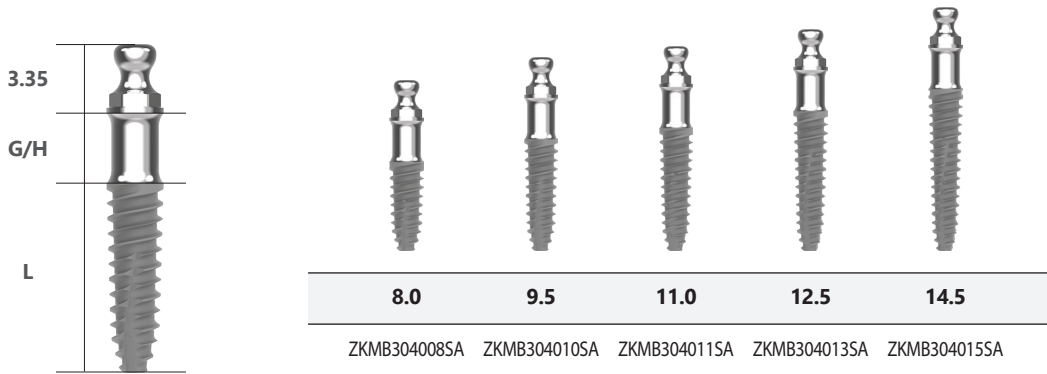

•Recommendable torque for Placement : Less than 40Ncm
 + FDA approval is under processing

Prosthetics Flow Diagrams for ZENEX Mini Post System

Cement Retained Restoration / Overdenture Retained Restoration

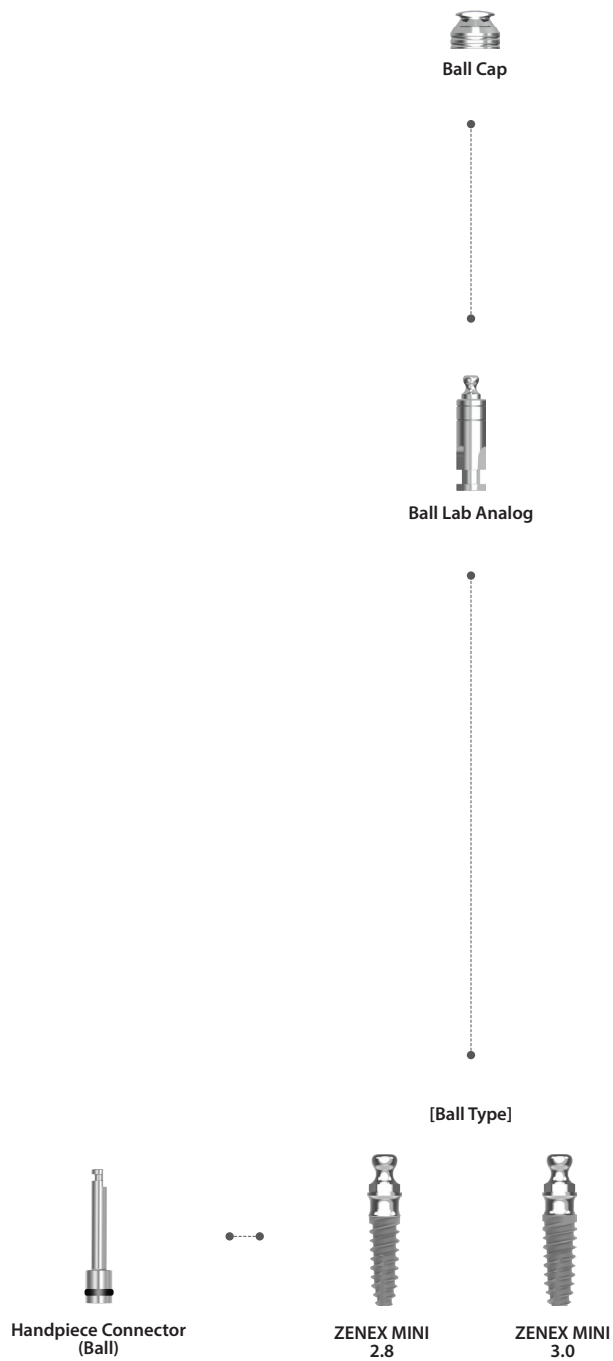

Ball Lab Analog



D	Ø4.0
-	ZATBLA

• The item which makes on the model of the actual Zenex Mini fixture (Ball Type).

Ball Cap O-Ring Set



D	Ø5.0
-	ZBAC

- 
Retentivity : 3~5Ncm (3ea of Black Ring included)
- Consist of O-ring
 - Allow maximum 20 degrees' path
 - For denture : Black color
 - Replace cycle : 6 months



HQ & Factory
1, 2Dong, 26-32, Suworam 4-gil, Seotan-myeon, Pyeongtaek-si,
Gyeonggi-do, Republic of Korea

Sales Office
#1013, #1014 Central-Tower II, 158, Magokseo-ro, Gangseo-gu, Seoul,
Republic of Korea
Tel. +82-2-3664-3666
Fax. +82-31-662-0658